

TECHNICAL DATA SHEET

BV-16SG White Surgical Gown

BV-16SG surgical gown, fully compliant with EN 14126:2003/AC:2004 standard, is made from non-woven PP+PE fabric providing AAMI PB70 level IV protection. Level IV gowns offer exceptional protection against fluids under high liquid pressure conditions. BV-16SG surgical gowns are designed to protect wearers from liquid splashes (Type PB[4]B)

KEY FEATURES

- Complies with EN 14126:2003/AC:2004
- AAMI PB70 level IV protection
- Single-use gowns
- Ties on the side eliminate the difficulty of tying behind the back
- Velcro closure on the neck provides perfect adjustment for full coverage and flexibility
- Cuffs made of white knitted material on the wrists area provide soft wrist closure for the comfort as well as ensuring a better and tighter fit
- Hot-sealed blue tapes provide extra protection against chemicals and liquids
- Anti-static treatment help reduce static build-up.
- Latex-free
- Non-Sterile
- Offers protection against biohazards and liquid splashes.

NON-STERILE

APPROVALS AND PROTECTION

EN 14126:
2003/AC:2004



EN 13034:
2005+A1:2009



TYPE 4
EN 14605:
2005+A1:2009



EN 1149-5:2018



EN ISO 13688:2013




Fulfill the requirements for Category III PPE. CE approved under Personal Protective Equipment Regulations, (EU) 2016/425



BELLEVITAE
www.bellevitae.com

SIZE

An appropriate size for the garment should be selected to allow sufficient movement for the wearer whilst maintaining a secure fit.

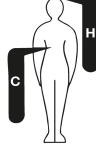

SIZE CHART

SIZE	HEIGHT	CHEST
S	116 cm	65 cm
M	118 cm	68 cm
L	120 cm	71 cm
XL	122 cm	74 cm
XXL	124 cm	77 cm
3XL	126 cm	80 cm

MATERIALS

SUIT	Non-woven PP+PE
WEIGHT	40 gr/m2

STORAGE AND DISPOSAL

- Store in dry, clean conditions in original packaging, away from direct sunlight, sources of high temperature, and solvent vapours.
- Replace garment if damaged, heavily contaminated or in accordance with local work practice or regulations.
- Expected shelf life is **five(5)** years from date of manufacture when stored as stated.
- Handle and dispose of contaminated garments with care and in accordance with applicable regulations.

WARNING AND LIMITATIONS

Make sure that the product is suitable for the application and fitted correctly. Product must never be altered or modified. Final determination as to the suitability of these products for a particular situation is the employer's responsibility. For full protection, hot-sealed tape along the seams has been applied. Electrostatic Dissipative (ESD) Protective clothing shall not be opened or removed whilst in presence of flammable or explosive atmospheres or while handling flammable or explosive substances. This information is subject to revision at any time.

PACKAGING.

Each gown is individually folded and packaged into a bag

LIMITED USE



Single use



Do not wash



Do not bleach



Do not iron



Do not dry clean



Do not tumble dry



Flammable

INTENDED USE

- For use in clinical and surgical operations.
- Intended to protect from the transfer of microorganisms and bodily fluids in moderate or high risk patient surgical situations.
- Typical applications may include contact with blood, bodily fluids and environmentally controlled cleanroom applications.